

Bee by Jenx

Styled specifically to appeal to very young children and their families, and developed in conjunction with therapists and carers, the Bee seating system by Jenx is an ideal playful and gentle introduction to supported seating. It provides an outstanding level of versatility through a comprehensive range of fixed and height adjustable bases, along with a wide variety of head, trunk, leg and foot support options. This means each Bee can be precisely tailored to the exact needs of each individual child. Adjustable back angle and tilt in space options, alongside the ability to accommodate babies as small as 4 months right through to children of around 5 years, mean the Bee really does tick every box for maximum functionality, flexibility and style!



INNOVATIVE FEATURE

Long Leg Sitting Board

Simply and firmly attached and removed with just two screws, the long leg sitting board can be used permanently or intermittently, both for very young children or those needing an extra stretch to the thighs or calves.

Angle-Adjustable Footplate

To accommodate or encourage dorsi or plantar flexion, adjustable quickly and simply with the turn of a button. Shown flipped up.

INNOVATIVE FEATURE

INNOVATIVE FEATURE

Back Angle Adjustment

From 10° recline through to 15° prone this feature can be used to create the required hip angle for the optimum seating position for each child.

INNOVATIVE FEATURE

Adjustable Seat Depth

To ensure the fit is exactly right for each child the seat depth is fully adjustable from 110-350mm.

INNOVATIVE FEATURE

Height Adjustment

4 different base options with a variety of fixed or adjustable heights allow the flexibility for each child to enjoy activities from floor level, to breakfast bar height.

Always Supportive

Every position a child adopts has an effect on their development and functional movement now and in the future.

Minimises fatigue and assists activity and functionality



Helps digestion



Encourages bone & muscle development



Assists neutral alignment



Proper support makes it easier for children to join in every day activities



Sitting well helps you breathe more easily



Support of pelvis, thighs & feet forms the foundation of stability



Provides opportunities to increase social skills



A stable base supports, develops and can provide greater independence and interaction with friends and family, resulting in more smiles for everyone!

INNOVATIVE FEATURE

Flip-away PU Support Pads

Included as standard for hip support and available in two sizes for optional trunk support. They have an easy, single-handed flip away action and are made from PU. These supports are fully flexible, easy to clean and impervious to fluids making them simple, hygienic and robust.

INNOVATIVE FEATURE

Adjustable Back Height

From 270-360mm as your child grows, the Bee grows with them!



Always Thoughtful

Outstanding features and benefits

Bright and Friendly

The diagnosis of a child with special needs is often very upsetting for families who may struggle to come to terms with all this can mean. Parents and carers can sometimes find it hard to accept the help and equipment their child needs at this difficult time. With this in mind the Jenx Bee has been carefully designed to offer the best in early intervention support in a gentle, bright and friendly way. Parents often comment on the huge difference the thoughtful design of the Jenx Bee has in helping them accept and use supported seating for their very young children.



Jenx Waistcoat Harness

Fantastic harnesses in a range of sizes to work across the entire range of Jenx seating products. Jenx put the same level of thoughtful design into their small products as their big ones, meaning that accessories are great pieces of kit in their own right as well as complementing the innovative design features of the main product. Providing firm and gentle trunk support to children with low tone, the Jenx waistcoat harness is easy for carers to use - made to attach, detach and flip away in seconds with a single hand. The zip is purposefully installed upside down, allowing access to abdominal feeding sites while maintaining good upper trunk control. A soft fabric surface against the child, with a wipe-clean front, make it both effective and practical for positioning support and for during feeding.

PU Support Pads

Polyurethane, or PU is a highly flexible and versatile material used across many elements in the Jenx range. It is impervious to fluids and therefore easy to clean and ideal from an infection control and re-issue point of view. Covers are available for many PU accessories and are recommended for additional comfort and to further extend the life of these items.



Height Adjustment

Young children typically spend their time playing on the floor, sitting in highchairs or working at low tables in their nursery or school environment. It can be a challenge to provide those children who need it the additional support in seating, whilst still enabling access to activities at every height necessary to facilitate engagement with their parents/carers and peers. The Bee range includes a choice of four different bases - from simple and affordable fixed-height options through to fully adjustable multi-height designs, there is something to suit every environment and budget, so no child misses out!

Jenx Multigrip™

The original Jenx Multigrip™ headrest is a unique and pioneering design, carefully developed to mimic the infinitely adjustable support offered by caring hands. Individual fingers on the Multigrip™ are adjusted to cushion, support and control almost any head position. Highly padded, removable and washable covers make Jenx Multigrip™ head supports a practical solution for almost any child. Multigrip™ head supports are available in a range of different length and height options, all including a choice of highly flexible fixing bracket, giving adjustment for height, depth and angle. Through this combination there is scope to block areas that may otherwise pose a hazard to children prone to getting their heads stuck beneath more conventional headrest designs.

Exclusively to the Jenx Bee, Multigrip™ technology is also available as a highly flexible shoulder protractor.

Visit jenx.com for additional product details.

Always Configurable

Accessories

The Jenx Bee's wide range of accessories and optional extras mean it can be highly customised for every child's needs. Visit jiraffe.org.uk for additional accessory details.

Bee Seat Unit

Suitable for children from as young as 4mths to around 5yrs, the Bee seat unit provides all the basics for a solid seating foundation, including:

- Padded seat and back
- Adjustable, flip away PU hip pads
- Adjustable footpost
- Angle-adjustable footplate
- Lapstrap with groin safety strap
- Push handle
- Flat headrest - this can also be used as an infill back pad for taller children.



Multigrip™ Headrests

The original and innovative Jenx Multigrip™ headrest is now available in 3 sizes for all Jenx seating systems. All Jenx Multigrip™ head supports have a choice of various highly flexible fixing brackets, all of which give adjustment for height, depth and angle.




Flat Head Rest/Upper Back Pad

Included as standard with every Bee seating unit and also available as an additional accessory, this is a simple headrest for children who require only very minimal head support. Can be used to extend the back pad height as the child and in addition to any other compatible Jenx headrest.



Multigrip™ Shoulder Protractor

A unique and innovative Jenx design, the trademarked Multigrip™ shoulder protractor is available exclusively on the bee. Cushioned fingers provide flexible support and it simply slots into place beneath the headrest, encouraging children to bring arms into midline.

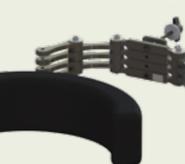

Slim Thoracic Support Pads

Available in a range of sizes with regular, or high/low rise brackets, these slimline trunk support pads are made from flexible and hygienic PU. Each size is available as a wrap around support with adjustable chest strap or shorter prompt support without strap. There is a single handed push button flip away action on all supports and they are independently adjustable for height and width. Can be offset to brace a scoliosis or accommodate other skeletal deformities.



Lumbar Pad

Simply Velcro into place to give additional lumbar support.



Additional Back Pad

A deeper 50mm backpad to reduce the seat depth for really small children. Simply Velcro into place behind the standard 25mm back pad - easily removed when no longer needed.



Long-Leg Sitting Extension

Easily and firmly attached and removed with just two screws, the long leg sitting board can be used permanently or intermittently for very young children or those needing an extra stretch to the thighs or calves. Adjustable support blocks firmly Velcro on. All other settings remain unchanged making the change from regular to long leg sitting quick and easy.



Foam Blocks

A pair of foam blocks which can be used to fill gaps or provide additional support or cushioning wherever it is needed. Covered in wipe clean PU fabric.



Waistcoat Harness

Provides additional trunk support to children with low tone. Attach, detach and flip away in seconds with a single hand – no tools required.



Armrests

Provided as a pair, the Jenx armrests are independently adjustable for height and angle. Clear markings make it easy to set the armrests at the correct height.



NB: Required for using a tray.

Nursery Tray

A useful surface for work, feeding or play, made from food grade plastic providing a hygienic and easy to clean surface. Contrasting tray edging gives a definitive boundary particularly useful for children with visual impairments and is also raised to allow toys or other items to be securely fastened in place. It's also available with reduced cut-out for babies and small children.



NB: Suction toys, anchors or grab bars will not work with this tray – to use these items with the Bee please see the nursery laminate tray.

Nursery Laminate Tray

A useful surface for work, feeding or play, made from composite laminate providing a hygienic and easy to clean surface. Contrasting tray edging gives a definitive boundary particularly useful for children with visual impairments and is also raised to allow toys or other items to be securely fastened in place. Perfect surface for use with suction toys, anchors or grab bars.



Padding on Tray Edge

Extra cushioning for the tray edge adjacent to the child. Attaches simply with Velcro and is covered in wipe-clean PU fabric.



Tray Elbow Blocks

Used to encourage arms into a midline position and/or to block excessive retraction.



Activity Frame

Great for encouraging head control and hand/eye co-ordination. Easy to attach and adjust. Comes with two clips for attaching toys or other items.



Grab Bar

Perfect for children who need to fix one hand to increase functionality.



Abduction Blocks

Available in two sizes providing different levels of abducted seating.



Knee Block

Easy to use knee block with fully adjustable knee cups which can be used to control legs and create abduction or adduction or control windsweeping.



X-Base

A compact, robust and cost effective base with high-low and tilt-in-space functionality. With lockable, swivel castors and a height range of 340-510mm. Tilt-In-Space 30° recline – 10° prone.



Y-Base

The most flexible base in the Jenx range, the Y-base offers a huge range of height adjustment alongside tilt-in-space functionality and lockable, swivel castors. Tilt-in-Space 30° recline – 10° prone.



Wooden High Base

Exceptional value, the Jenx Wooden High Base for the Bee Seating System offers fixed height seating (550mm) which is compact, economical and stylish. It features lockable swivel castors and Tilt-In-Space with a recline to 20°.



Wooden Low Base

Exceptional value, the Jenx Wooden Low Base for the Bee Seating System offers fixed height seating (400mm) which is compact, economical and stylish. It features lockable swivel castors and Tilt-In-Space with a recline to 20°.



Footplate Pad

Attach to the footplate to give a softer surface for children who don't need more structured foot positioning.



Sandals

Available in 4 different sizes to accommodate any size of foot, even when wearing splints, Jenx sandals are interchangeable across the entire range of Jenx seating and standing products.



Sandal Raising Blocks

Available for every size of Jenx sandals, the raising block kits give the sandals that extra bit of height where needed! Use to accommodate for leg length discrepancy or use to raise the height of the footplate on any product to maximise the available growth.



NB: One kit needed per individual sandal.

Booties

Available for the smallest Jenx seating systems, booties provide gentle and comfortable support for tiny feet even without shoes. They are soft and extra padded with no fiddly straps, press studs or buckles.



Always Adaptable



Jenx Bee Technical Specifications

	Units	Bee NM01
Age Range (approx)	mths/yrs	4mths - 5yrs
Max User Weight	kg/lbs	25 / 55
Back Angle (recline/prone)	degrees	10° - 15°
Back Pad Height (seat to top min/max)	mm/inch	270 - 360 / 10½ - 14
X-Base (footprint)	mm/inch	610 - 620 / 24 - 24½
Y-Base (footprint)	mm/inch	650 - 780 / 25½ - 30¾
Wooden Base (footprint)	mm/inch	500 - 600 / 19¾ - 23¾
Chest Width (min/max)	mm/inch	150 - 270 / 6 - 10½
Hip Width (min/max)	mm/inch	130 - 250 / 5 - 9¾
Seat Height From Floor (X-Base min/max)	mm/inch	340 - 510 / 13¼ - 20
Seat Height From Floor (Y-Base min/max)	mm/inch	285 - 750 / 11¼ - 29¾
Seat Height From Floor (Wooden Low Base/High Base)	mm/inch	400 - 550 / 15¾ - 21¾
Seat Depth	mm/inch	110 - 350 / 4¼ - 13¾
Lower Leg Length	mm/inch	95 - 345 / 3¾ - 13½
Depth With Long Leg Conversion (min/max)	mm/inch	370 - 560 / 14½ - 22
Tilt-In-Space Angle (X and Y Base recline/prone)	degrees	30° - 10°
Tilt-In-Space Angle (Wooden Base Recline)	degrees	20°

Always here

To arrange a demonstration or assessment with the Jenx Bee seating system, just get in touch!

My name is:

Please ring me on:

Or email me at:

Seating

Standing

Sleeping

Therapy

JENX[®]
Positioning for life

Bee by Jenx

A fully adjustable, modular seating system specifically designed for young children from 4 months – 5 years



Always Supportive...

Jenx Limited, Wardsend Road, Sheffield S6 1RQ, United Kingdom
tel: +44 (0) 114 285 3376
email: export@jenx.com
web: www.jenx.com

JENX[®]
Positioning for life

Jenx Limited, Wardsend Road, Sheffield S6 1RQ, United Kingdom
tel: +44 (0) 114 285 3376
email: export@jenx.com
www.jenx.com

JENX[®]
Positioning for life